

<p>OMRON HEALTHCARE EUROPE B.V. Scorpions 35, 2132 LF Hoofddorp, THE NETHERLANDS</p>	<p>Reverendur www.omron-healthcare.com</p>
<p>OMRON HEALTHCARE UK LTD Opal Drive, Fox Hills, Milton Keynes, MK15 9GQ, UK</p>	<p>Sustisaires www.omron-healthcare.com/distributors</p>
<p>OMRON SANTE FRANCE SAS www.omron-healthcare.com/distributors</p>	<p>www.omron-healthcare.com/distributors</p>
<p>OMRON MEDIZINTECHNIK HANDELSGESELLSCHAFT MBH www.omron-healthcare.com/distributors</p>	<p>www.omron-healthcare.com/distributors</p>
<p>OMRON HEALTHCARE S.r.l. Via Marziale Central, 90/Fc 25017 Lonato del Garda (BS) Italy</p>	<p>Fabritant www.omron-healthcare.com</p>


Issue Date / Date de publication: 2022-11-16
IM1-NE-C301-E-02-11/2022
3A44518 rev.01

Instruction Manual

DuoBaby (NE-C301-E)


Compressor Nebulizer with Nasal Aspirator

Nébuliseur à compresseur avec aspirateur nasal

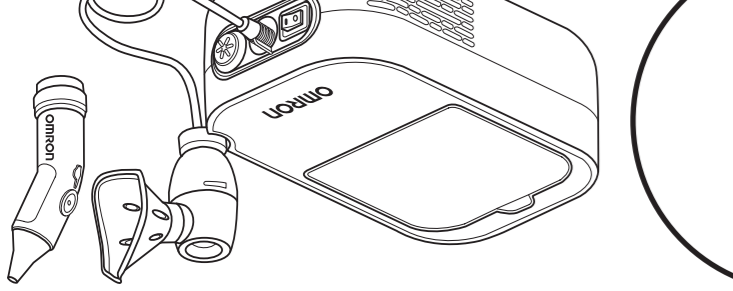


OMRON


Know your unit

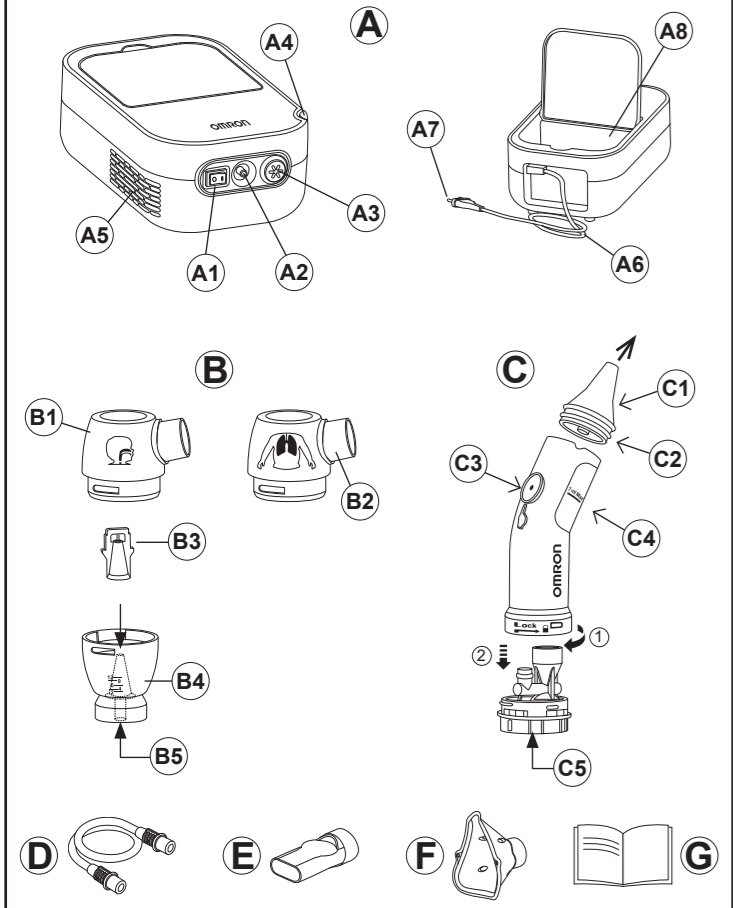


How to use the nebulizer kit

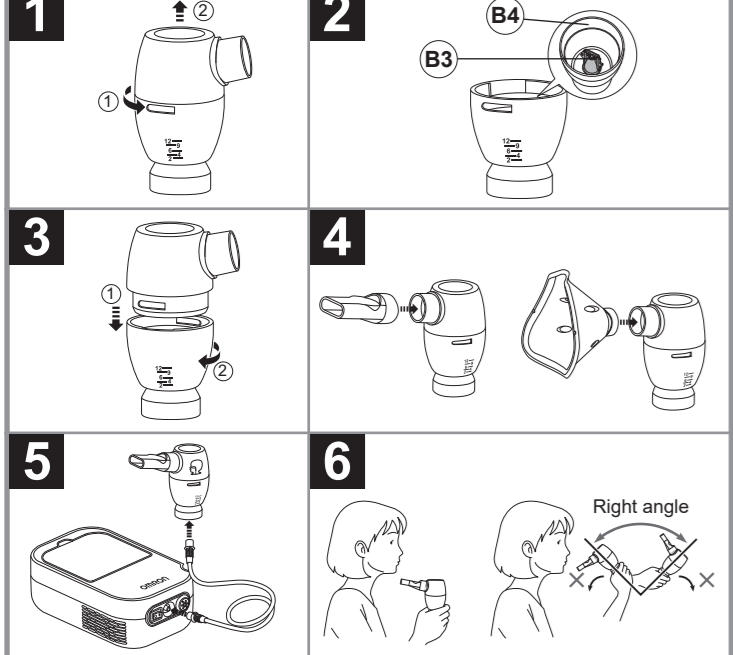

How to use the Nasal Aspirator

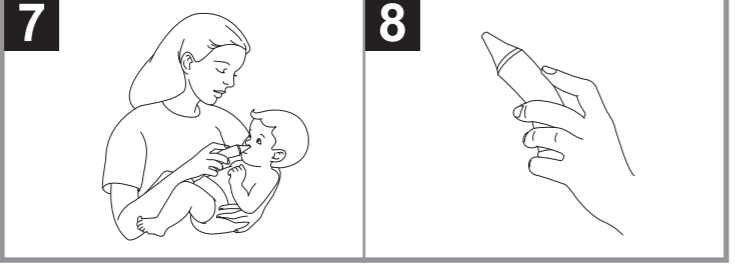

Replacing the Air Filter

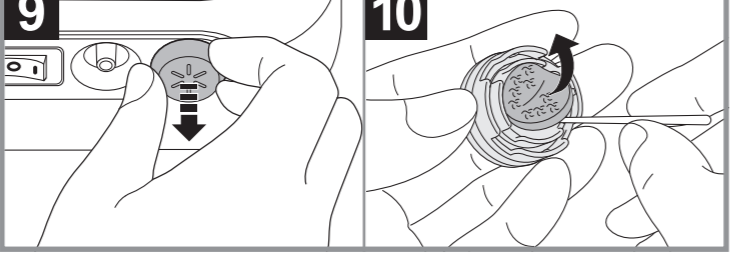

Intended use

Medical Purpose
The intended purpose of this product is:
- To inhale medication for respiratory disorders.
- To remove mucus from the nasal cavities to ease breathing.

The compressor is intended for multi-patient reuse, while the nasal aspirator, the nebulizer kit and its attachments are single-patient reuse.

Intended User
• Patient and/or caregiver under the direction of legally certified medical experts, such as doctor, nurse and therapist or qualified medical experts for home treatment.
• The user should also be capable of understanding general operation of DuoBaby and the content of this instruction manual.

Intended Patients
Persons from the age of 1 month who suffer from upper airway conditions (such as blocked nose, cold, allergies) and/or lower respiratory conditions (such as asthma, bronchitis, bronchiolitis).

Environment
This product is intended for use in a general household.

Durable period
Durable periods are as follows, provided the product is used to nebulize 3 times a day for 10 minutes each time at room temperature (23°C). Durable period may vary depending on usage and environment.

Compressor (Main unit): 5 years
Nebulizer Kit - Mouthpiece - Air Tube - Mask - Nasal Aspirator: 1 year
Air Filter: 70 applications

Precautions for use
Warnings and cautions described in the instruction manual should be observed.


How the nebulizer kit works

The medication that is pumped up through the medication channel is mixed with compressed air which is generated by a compressor pump. The compressed air mixed with medication is turned into fine particles and sprayed when in contact with the baffle on the vaporiser head.

How the nasal kit works

The compressed air generated by the compressor pump, passing through a narrowing, thanks to the Venturi effect, creates a vacuum which sucks the mucus into the collection chamber.

Know your unit

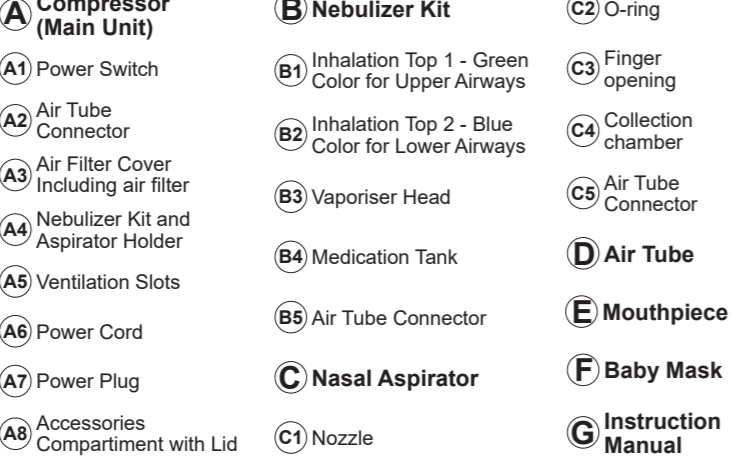

Important safety instructions

- Warning:** Indicates a potentially hazardous situation which, if not avoided, could result in serious injury.
- The nasal aspirator, the nebulizer kit and its attachments are intended for use by a single patient. Do not use by multiple patients, use separate accessories for each patient.
 - Clean the parts after each use and disinfect them once a week during normal use.
 - Do not use the device for other use than nebulization for human aerosol therapy and for removing the mucus.
 - Manufacturer will not be held responsible for inappropriate use of the device.
 - Use only original OMRON accessories as indicated in this instruction manual.
 - The correct functioning of the device may be affected by electromagnetic interference caused by malfunctioning television sets, radios, etc. If this occurs, move the source of interference away from the device or use a different plug.
 - For type, dose, and regime of medication follow the instructions of your doctor or respiratory therapist. Before use check that you are using the correct Inhalation top.
 - If you feel or see anything unusual during use, stop using the device immediately and consult your doctor.
 - Do not use only water in the nebulizer for inhaling purposes.
 - Keep the device out of the reach of unsupervised infants and children. The device may contain small pieces that can be swallowed.
 - In view of their length, the power cord and air tube could constitute a strangulation hazard.
 - Do not store the air tube while there is moisture remaining inside.
 - Always dispose of any remaining medication after use, and use fresh medication each time.
 - Do not use in anaesthetic or ventilator breathing circuits.
 - Never block the ventilation slots. Never place the device where the ventilation slots may be obstructed during operation.
 - Never leave the nebulizer kit with medication near the ventilation slots, this is to avoid medication entering the device.
 - Store the parts in a clean location for avoiding infection.

- (Risk of Aspiration)**
- Never plug or unplug the device with wet hands.
 - The compressor and power cord are not waterproof. Do not spill water or other liquids on these parts. If liquid does spill on these parts, immediately unplug the power cord and wipe off the liquid.
 - Do not immerse the compressor in water or other liquid.
 - Do not use or store the device in humid locations.
 - Do not operate the device with a damaged power cord or plug.
 - Keep the power cord away from heated surfaces.
 - Always unplug the device from the power outlet after use and before cleaning.
 - Plug the device into the appropriate voltage outlet. Do not overload power outlets or use extension cords.
 - Do not misuse the power cord.
 - Do no wind the power cord around the compressor.
 - Be sure to use the compressor at a place where the power plug is easily accessible during treatment.

- (Cleaning and disinfecting)**
Observe the rules below when cleaning or disinfecting parts. Failure to observe these rules may result in damage, inefficient nebulization or infection. For the instructions, refer to "Cleaning and disinfecting" section.
- Clean and disinfect the nebulizer kit, baby mask, mouthpiece and nasal aspirator before using them when as follows:
 - the first time after purchase
 - if the device has not been used for a long period of time
 - Clean the parts after each use and disinfect them once a week during normal use. The parts are expected to withstand approx. 50 disinfection cycles.
 - Be sure to wash or wipe the parts after use, and ensure that they are thoroughly dried and stored in a clean location.
 - Do not leave the cleaning solution in the parts.
 - Do not use a microwave oven, dish dryer or hair dryer to dry the device or the parts.
 - Do not use an autoclave, EOG gas sterilizer or low temperature plasma sterilizer.
 - When disinfecting parts by boiling (for at least 5 minutes), make sure that the container does not boil dry. Otherwise it may also result in fire.
- Caution:** Indicates a potentially hazardous situation which if not avoided, may result in minor or moderate injury, or physical damage.

- General Safety Precautions:**
- Inspect the device and parts before using them each time, and check that there are no problems. In particular, be sure to check the following:
 - That the nozzle of either the nebulizer kit or aspirator or air tube are not damaged or blocked.
 - That the compressor operates normally.
 - When using this device, there will be some noise and vibration caused by the device. This is normal and does not indicate a malfunction.
 - Do not use if the device or any parts are damaged. Replace damaged parts.
 - Do not subject the device or the parts to any strong shocks such as dropping the device on the floor.
 - Do not insert objects inside the compressor.
 - Do not disassembly or attempt to repair the compressor or power cord.
 - Do not leave the device or its parts in extreme hot or cold temperature, or under direct sunlight.

Keep these instructions for future reference.

- (Nebulizer: Usage)**
- Provide close supervision when this device is used by, on, or near children, invalids or persons who cannot express their consent.
 - Make sure that the parts are attached or installed correctly.
 - Make sure that the air filter is clean and correctly attached. If the air filter has been used for more than 70 applications, replace it with a new one.
 - Do not tilt the nebulizer kit at an angle of greater than 60 degrees in all directions or shake it while in use.
 - Do not use or store the device while the air tube is creased.
 - Do not add more than 12 ml of medication to the medication tank.
 - Do not carry or leave the nebulizer kit while the medication tank contains medication.
 - This product should not be used on patients, who are unconscious, sleeping, are drowsy or are not breathing spontaneously.
 - Do not distort the vaporiser head, or jab the nozzle of the medication tank with a pin or any sharp object.
 - The compressor may become warm during operation, therefore do not touch it other than turning it off after nebulizing.
 - To avoid medication residue on the face, be sure to clean the patient's face after removing the mask.

- (Nasal Aspirator: Usage)**
- The nasal aspirator should be used exclusively for removing the nasal mucus from a baby's nose. It cannot be used for any other purpose.
 - Do not obstruct the holes on the bottom of the nasal aspirator.
 - If you notice rashes or redness appearing, inside or around the nose of the baby during or after using the device or anything unusual during use, stop using the device and consult a physician.
 - Only use the original parts as included in the product.
 - Do not attach the nasal aspirator to any other device.

User Guidance about the Nebulizer Kit

The DuoBaby nebulizer is a medical device that generates an aerosol with adjustable characteristics to adapt to the patient's respiratory pathology

Inhalation Top 1 – Green Color for Upper Airways:	Treats conditions such as:	Inhalation Top 2 – Blue Color for Lower Airways:	Treats conditions such as:
rhinitis, sinusitis, pharyngitis, tonsillitis, and laryngitis.	asthma, bronchitis, bronchiolitis, bronchiectasis and bronchopneumonia.	asthma, bronchitis, bronchiolitis, bronchiectasis and bronchopneumonia.	asthma, bronchitis, bronchiolitis, bronchiectasis and bronchopneumonia.
Particle size (MMAD) is approximately 9.3 µm.	Particle size (MMAD) is approximately 3.9 µm.	Particle size (MMAD) is approximately 3.9 µm.	Particle size (MMAD) is approximately 3.9 µm.

How to use the nebulizer kit

- Make sure that the power switch is in the off (O) position and the power plug is disconnected from the power outlet.
- Remove the inhalation top from the medication tank **1**.
 - Rotate the inhalation top anti-clockwise
 - Lift the inhalation top out from the medication tank.
- Verify the presence of the vaporiser head (B3) inside the medication tank.
- Add the correct amount of prescribed medication into the medication tank **2**.
- Choose the appropriate inhalation top by reading the "User Guidance". To put the inhalation top back onto the medication tank **3**.
 - Lower the inhalation top onto the medication tank.
 - Rotate the inhalation top clockwise.
- Attach the Baby mask or mouthpiece to the nebulizer kit tightly **4**.
- Attach the air tube. While twisting the air tube plug slightly, push it firmly into the air tube connector **5**.
- Hold the nebulizer kit as indicated in figure **6**. Follow the instructions of your doctor or respiratory therapist.

Caution: Do not tilt the nebulizer kit at an angle of greater than 60 degrees in all directions. Medication may flow into the mouth or it may result in ineffective nebulization.
- Plug the power plug into a power outlet.

- Turn the power switch to the on (I) position. As the compressor starts, nebulization begins and aerosol is generated. Inhale the medication. Exhale through the nebulizer kit.

- When treatment is completed, turn the power off and unplug the compressor from the power outlet.

- Clean and disinfect all (disassembled) parts according to the "Cleaning and disinfecting" section.

How to use the Nasal Aspirator

Dry and very thick mucus inside the nose may be dissolved by saline solutions. After adding saline solution into each nostril of the baby's nose to help loosen nasal secretions, wait a few minutes and then remove the mucus running out with a tissue. Only after removing the running mucus outside the nostrils the nasal aspirator may be used.

- Make sure that the power switch is in the off (O) position and the power plug is disconnected from the power outlet.
- Assemble the nasal aspirator. Make sure that the O-ring is properly attached. Refer to the drawing "Nasal Aspirator".
- Connect the air tube to the nasal aspirator and to the main unit, while twisting the air tube plugs slightly, push them firmly.
- Keep your baby in a comfortable position and make sure that their head is lifted as shown on the figure **7**.
- Plug the power plug into a power outlet.
- Turn the power switch to the on (I) position. Place your index finger over the opening on the nasal aspirator, as shown on the figure **8** to begin aspiration.

Caution: Before using the nasal aspirator, try it on the hand to see if it sucks. When functioning properly, there must be a slight suction from the nozzle. If air is exiting from the nozzle, do not continue using the aspirator. While using, the aspirator needs to be in the upright position, with the finger opening on top (see figure by above). This prevents the sucked mucus from entering parts other than the collection chamber.

- Place the nozzle at the entrance of one of the nostrils of the baby. (do not insert inside the nostril). Use the aspirator for 2 seconds for 4 times per nostril.
- When treatment is completed, turn the power off and unplug the compressor from the power outlet.
- Disconnect the aspirator from the air tube and disassemble it into 3 parts as shown in "Know your unit"; clean or disinfect according to the section below.

Cleaning and disinfecting

Cleaning
Clean the parts after each use to remove residual medication or mucus. This will prevent inefficient nebulization or aspiration and reduce the risk of infection. Do not use a brush or pin to clean.

- Washable parts: **Nebulizer kit, Baby Mask (PVC), Nasal Aspirator, Mouthpiece**
Wash them in warm water and mild, neutral detergent. Rinse them thoroughly with clean hot tap water, gently tap to remove excess water and allow to air dry in a clean place.

- Non-washable parts: **Compressor, Air Tube (PVC)**
Firstly, make sure that the power plug is unplugged from the power outlet. Wipe clean with a soft cloth, moistened with water or mild, neutral detergent.

Air Filter
Do not wash or clean the air filter. Refer to "Changing the air filter" section.

Disinfecting
Always disinfect the parts before using the product for the first time and after the product has not been used for an extended period of time. Disinfect the parts once a week. The parts are expected to withstand approx. 50 disinfection cycles. If the parts are heavily stained, replace them with new ones. To select a method for disinfection, refer to the table below.

Disinfecting method	Parts	Mouthpiece	Nasal Aspirator	Nebulizer kit	Baby Mask*	Air Tube (Phthalate free)	Air Filter
A See notes below	Disinfecting ethanol					X	X
	Sodium hypochlorite					X	X
	Milton*					X	X
	Quaternary ammonium					X	X
	Osvan*					X	X
B	Chlorhexidine					X	X
	Hibitane*					X	X
	Amphoteric Surfactant					X	X
B	Tego*					X	X
	Boiling				X	X	X

- Notes:**
- Use a commercially available disinfectant. Follow the instructions provided by the disinfectant manufacturer.
 - Never clean with benzene, thinner or a flammable chemical.
 - After having washed the Nebulizer Kit and/or Nasal Aspirator with the available disinfectant, please rinse with clean warm water and dry thoroughly before using.
 - Applicable parts may be boiled for at least 5 minutes.
 - After boiling, carefully remove the parts, shake off excess water and allow to air dry in a clean environment.

Removing condensation from the air tube

- If there is moisture remaining in the air tube, be sure to follow the procedure below.
- Make sure that the air tube is connected to the air connector on the compressor.
 - Remove the air tube from the nebulizer kit or the nasal aspirator.
 - Turn on the compressor and pump air through the air tube to expel the moisture.

Changing the air filter

If the air filter has changed colour, or has been used for more than 70 applications, replace it with a new one. (See Other Optional / Replacement Parts).

- Remove the air filter cover as shown in figure **9**.
- Remove the old air filter with a sharp object, such as a toothpick, and insert the new air filter **10**.
- Put the air filter cover back in place.

- Notes:**
- Use only OMRON air filters designed for this device. Do not operate without a filter.
 - Do not attempt to wash or clean the filter. If the air filter becomes wet, replace it. Damp air filters can cause blockages and may increase risk of infection.
 - There is no front or back orientation for the air filters.
 - Check that air filters are clean and free of dust before inserting them.

Troubleshooting for: main unit, nebulizer kit and nasal aspirator

Check the following if the device should fail before or during operation.

Problem	Cause	Remedy
Nothing happens when the power switch is turned.	Is the power plug plugged into a power outlet?	Check that the plug is inserted in a power outlet. Unplug then reinsert the plug if necessary.
No nebulization/ aspiration or low nebulization, when the power is on.	Is there sufficient medication in the medication tank? (only nebulizer kit)	Add the correct amount of medication to the medication tank.
	Is the vaporizer head missing or not assembled correctly? (only nebulizer kit)	Attach the vaporizer head correctly.
Is the nebulizer kit or nasal aspirator assembled correctly?	Is the nebulizer kit or nasal aspirator assembled correctly?	Assemble the nebulizer kit or nasal aspirator correctly.
	Is the nozzle for the nebulizer kit or aspirator blocked?	Make sure that the nozzle is free of blockages, do not pinch the nozzle.
Is the nebulizer kit or nasal aspirator tilted at a sharp angle?	Is the nebulizer kit or nasal aspirator tilted at a sharp angle?	Make sure that the nebulizer kit is tilted at an angle less than 60 degrees, or the aspirator is upright.
	Is the air tube connected correctly?	Make sure that the air tube is correctly connected to the compressor and nebulizer kit or nasal aspirator.
Is the air tube folded, blocked or damaged?	Is the air tube folded, blocked or damaged?	Make sure that the air tube does not contain kinks or is not damaged or blocked.
	Is the air filter dirty?	Replace the air filter with a new one.
The compressor is abnormally loud.	Is the air filter cover attached correctly?	Attach the air filter cover correctly.
The compressor is very hot.	The compressor is insufficiently ventilated.	Do not cover the compressor during use.
		Do not block the ventilation slots.
The compressor suddenly stops during use.	Do not use the device in an environment with temperatures higher than 40°C.	
	The thermal cut-out has shut the device down due to overheating.	See note below.
Note:		
Do not attempt to repair the device. Do not open and/or tamper with the device. No parts of the device are user serviceable. Return the device to an authorized OMRON retail outlet or distributor.		



Please report to the manufacturer and the competent authority of the Member State in which you are established about any serious incident that has occurred in relation to this device.

Technical data

Product Description:	Compressor Nebulizer with Nasal Aspirator
Model (code):	DuoBaby (NE-C301-E)
Rating:	230 V ~ 50 Hz
Power Consumption:	140 VA (70W)
Operating Temperature/ Humidity/Air Pressure:	+10°C to +40°C / 10% to 95%RH / 690 - 1060 hPa
Storage and Transport Temperature/Humidity/ Air Pressure:	-20°C to +70°C / 5% to 95%RH / 690 - 1060 hPa
Classifications:	Type BF (Applied part): Mouthpiece, Mask and Nozzle of Nasal Aspirator Class II (Protection against electric shock) IP21 (Ingress Protection)
IP Classification:	IP21
Weight:	Approx. 1.7 kg (compressor only)
Dimensions:	Approx. 170 (W) x 280 (D) x 105.5 (H) mm
Contents:	Compressor, Nebulizer Kit, Air tube, Mouthpiece, Baby Mask, Nasal Aspirator, Instruction Manual.

- = Class II equipment
- = Type BF applied part
- = IP21 = Ingress Protection
- = Power off
- = Serial number
- = Consult the instructions for use
- = ~ = Alternating current
- = Power on
- = Medical Device
- = Consult the instructions for use
- = Date of Manufacture
- = UDI = Unique device identifier

	Nebulizer Upper Airways	Nebulizer Lower Airways
Nebulization Rate (by weight loss):	Approx. 0.4 ml/min (NaCl 0.9%)	Approx. 0.2 ml/min (NaCl 0.9%)
Appropriate Medication Quantities:	1 ml minimum – 12 ml maximum	
Sound:	Noise level (at 1m distance): Approx. 63dB	
Nebulization Rate (by weight loss):	Approx. 0.4 ml/min (NaCl 0.9%)	Approx. 0.2 ml/min (NaCl 0.9%)
Aerosol Output Rate**:	Approx. 0.14 ml/min	Approx. 0.07 ml/min
Percentage of fill volume emitted in 1 minute**:	Approx. 4.7%	Approx. 2%
Aerosol Output**:	Approx. 0.47 ml	Approx. 0.59 ml
Residual volume**:	Approx. 1.6 ml	Approx. 1.2 ml
Particle Size (MMAD)**:	Approx. 9.3 µm	Approx. 3.9 µm
GSD**:	Approx. 2.8	Approx. 2.3
Respirable fraction < 5 µm**:	Approx. 23%	Approx. 61%
Aerosol fraction < 2 µm**:	Approx. 6%	Approx. 22%
Aerosol fraction from 2 µm to 5 µm**:	Approx. 17%	Approx. 39%
Aerosol fraction > 5 µm**:	Approx. 77%	Approx. 39%
MMAD = Mass Median Aerodynamic Diameter GSD = Geometric standard deviation		
** 3ml, 0.1% Salbutamol in 0.9% NaCl solution.		

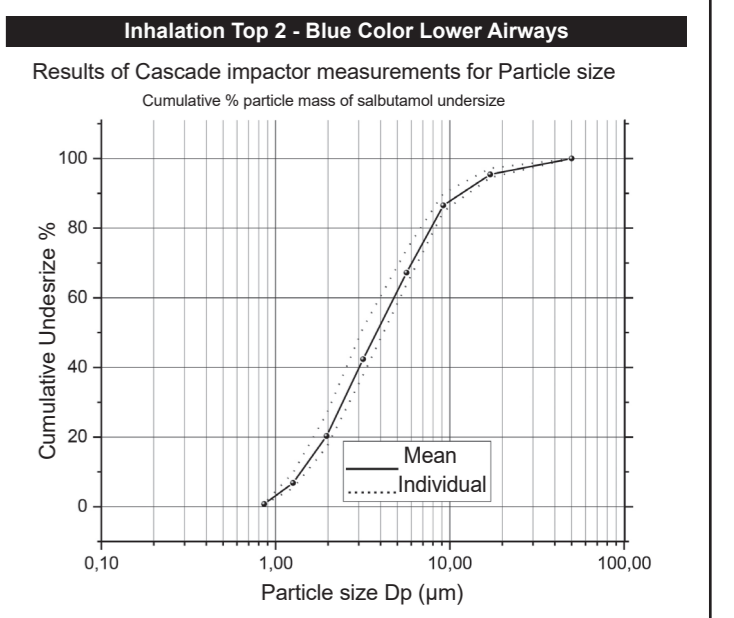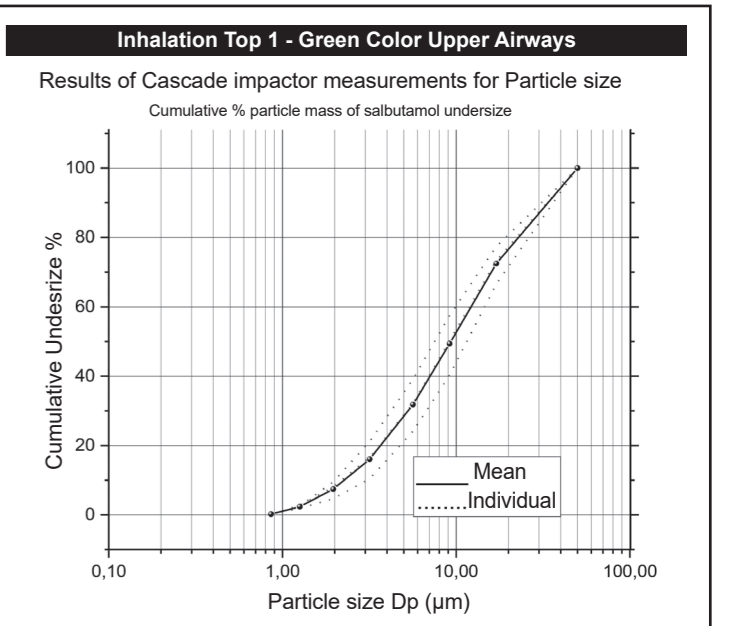

Nasal Aspirator
Suction Pressure: Minimum -0.15 bar - Maximum -0.2 bar
Air Flow: Approx. 4 l/min
Maximum fill volume of collection chamber: 7 ml

- General Notes:**
- Subject to technical modification without prior notice.
 - The device may not work if the temperature and voltage conditions are different to those defined in the specifications.
 - This device (compressor, nebulizer kit, air tube, mouthpiece and Baby mask) complies with the requirements of the European Standard EN ISO 27427:2019 Anaesthetic and respiratory equipment - Nebulizing systems and components.
 - The nebulizer performance may vary with drugs such as suspensions or high viscosity solutions. See drug supplier's data sheet for further details.
 - The particle size distribution and aerosol output may vary by combination of product, medication and ambient conditions such as temperature, humidity and atmospheric pressure.
 - The disclosures for nebulizer performance are based upon testing that utilizes adult ventilatory patterns and are likely to be different from those stated for paediatric or infant populations.
 - IP classification is degrees of protection provided by enclosures in accordance with IEC 60529. The compressor is protected against solid foreign objects of 12.5 mm diameter and greater such as a finger, and against vertically falling water drops which may cause issues during a normal operation.
 - See website of OMRON HEALTHCARE EUROPE for complete up to date information. URL: www.omron-healthcare.com

Warranty

Thank you for buying an OMRON product. This product is constructed of high quality materials and great care has been taken in its manufacturing. It is designed to give you a high level of comfort, provided that it is properly operated and maintained as described in the instruction manual. This product is guaranteed by OMRON for a period of 3 years after the date of purchase. The proper construction, workmanship and materials of this product is guaranteed by OMRON. During this period of guarantee OMRON will, without charge for labour or parts, repair or replace the defect product or any defective parts.

The guarantee does not cover any of the following:

- Transport costs and risks of transport
- Costs for repairs and / or defects resulting from repairs done by unauthorised persons.
- Periodic check-ups and maintenance
- Failure or wear of optional parts or other attachments other than the main device itself, unless explicitly guaranteed above.
- Costs arising due to non-acceptance of a claim (those will be charged for).
- Damages of any kind including personal caused accidentally or from misuse. Should guarantee service be required please apply to the dealer whom the product was purchased from or an authorised OMRON distributor. For the address refer to the product packaging / literature or to your specialised retailer.

If you have difficulties in finding OMRON customer services, contact us for information. Repair or replacement under the guarantee does not give rise to any extension or renewal of the guarantee period. The guarantee will be granted only if the complete product is returned together with the original invoice / cash ticket issued to the consumer by the retailer. OMRON reserves the right to refuse the guarantee service if any unclear information has been given.

Optional Medical Accessories

Product description	Model
DuoBaby Nebulizer accessory set (contents: nebulizer kit (includes two inhalation tops), baby mask, adult mask, mouthpiece, 3 air filters, air tube)	NEB2019
DuoBaby Aspirator accessory set, (contents: nasal aspirator, air tube)	NEB7000

Other Optional/Replacement Parts

Product description	Model
Air Filter set (contents: 3 pieces)	3AC408

DISPOSAL PROCEDURE (Dir. 2012/19/EU-WEEE)
This product is not to be treated as regular household waste but must be returned to a collection point for recycling electric and electronic devices. Further information is available from your municipality, your municipality's waste disposal services, or the retailer where you purchased your product.

Important information regarding Electro Magnetic Compatibility (EMC)
NE-C301-E manufactured by 3A HEALTH CARE S.r.l

